


MINDRAY

BC-3200 Auto Hematology Analyzer



Closed Tube Sampling Technology

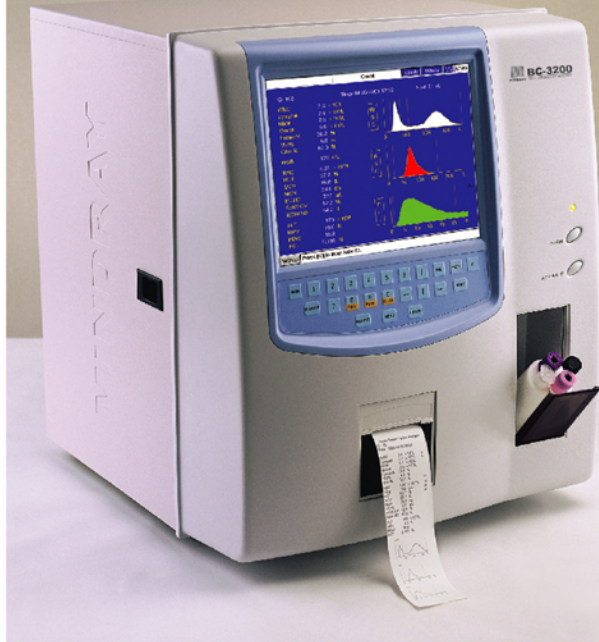

www.mindray.com

BC-3200 Auto Hematology Analyzer



- Closed tube sampling
- Sampling with 4 position tube holder
- Micro sampling of 13μl
- 3-part differential, 19 parameters +3 histograms
- Up to 60 samples per hour
- Large color LCD display
- Data storage of 35,000 samples with histogram


Flexible and Easy To Use

- Fully automatic analysis cycle
- Intelligent alarm to avoid wrong tube
- Automatic internal and external needle cleaning



High Reliability and Safety

- Sampling on closed tube to avoid biohazard
- Automatic protection when chamber door is open
- Cyanide-free reagent to protect environment

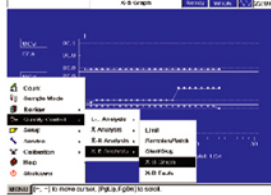

Total Solution for QC and Calibration

- Comprehensive QC programs including L-J, $\bar{X}$ -R and X-B analysis
- Up to 9 QC file - 31 runs/file
- Fully automatic calibration
- Unique calibration program with fresh blood



Powerful Customizable Software Design

- Huge storage of 35,000 patient results with histogram
- Convenient patient data review
- Programmable multi-format printout
- Fully automatic daily maintenance program

BC-3200 Auto Hematology Analyzer

Technical Specifications

Parameters

WBC, Lymph#, Mid#, Gran#, Lymph%, Mid%, Gran%, RBC, HGB, HCT, MCV, MCH, MCHC, RDW-CV, RDW-SD, PLT, MPV, PDW, PCT and Histogram for WBC, RBC, PLT

Principles

Electrical resistance for counting and SFT method for hemoglobin

Performance

Parameter	Linearity Range	Precision(CV%)
WBC($10^9/L$)	0.3-99.9	2.5 (7.0-15.0)
RBC($10^{12}/L$)	0.20-8.00	2 (3.5-6.0)
HGB(g/L)	10-250	1.5 (110-180)
MCV(fL)	0-5	0.5 (80.0-110.0)
PLT($10^9/L$)	10-999	5 (150-500)

Sample Volume

Prediluted	20μl
Whole Blood	13μl

Sample processing

On closed tube with 4 position tube holder (QC, 1.5mL, 3mL, 5 mL)

Aperture Diameter

WBC	100μm
RBC/PLT	70μm

Dilution Ratios

	Whole blood	Capillary Blood
WBC	1:308	1:546
RBC/PLT	1:44872	1:45827

Throughput

60 samples per hour

Display

Large Color LCD display
Resolution: 800x600

Alerts

45 sets of alert message

Menu

Count, Sample Mode, Review, Quality Control, Setup, Service, Calibration, Help, Shutdown

Carryover

WBC, RBC, HGB < 0.5%, PLT < 1%

Input/Output

RS232x2, 1 Parallel printer(optional), 1 bar code scanner(optional), 1 keyboard interface

Printout

Thermal recorder, 50mm wide paper, various printout formats, printer optional

Operating Environment

Temperature: 15°C~30°C
Humidity: 30%-85%

Power Requirement

100-240V~50/60Hz

Dimension

457mm(H)x 391mm(W)x 415mm(D)

Weight

25Kg



MINDRAY™

Headquarters: Mindray Building, Keji 12th Road South, High-tech Industrial Park, Nanshan, Shenzhen 518057, P. R. China

Tel: +86 755 28582492, 28582888 Fax: +86 755 28582500, 28582501
E-mail: info-market@mindray.com Web site: www.mindray.com



DISTRIBUTOR: